

Proclaim™

Quick Start Guide



For how-to videos, visit
[proclaimhealth.com/
instructions](https://proclaimhealth.com/instructions) or scan
the code.




TIPS FOR A GREAT EXPERIENCE:

- ♦ Start on pressure level 1 and increase to **level 5** over time
- ♦ Use with **warm water**
- ♦ Add Proclaim Oral Infusion

Welcome To A Refreshing, One-Of-A-Kind Experience

As the new owner of a Proclaim Custom-Jet Oral Health System, you can now elevate your oral care at home in just 7 seconds a day. Learn about your Proclaim system below, and then follow the simple steps to get started.


Get To Know Your Proclaim Custom-Jet Oral Health System

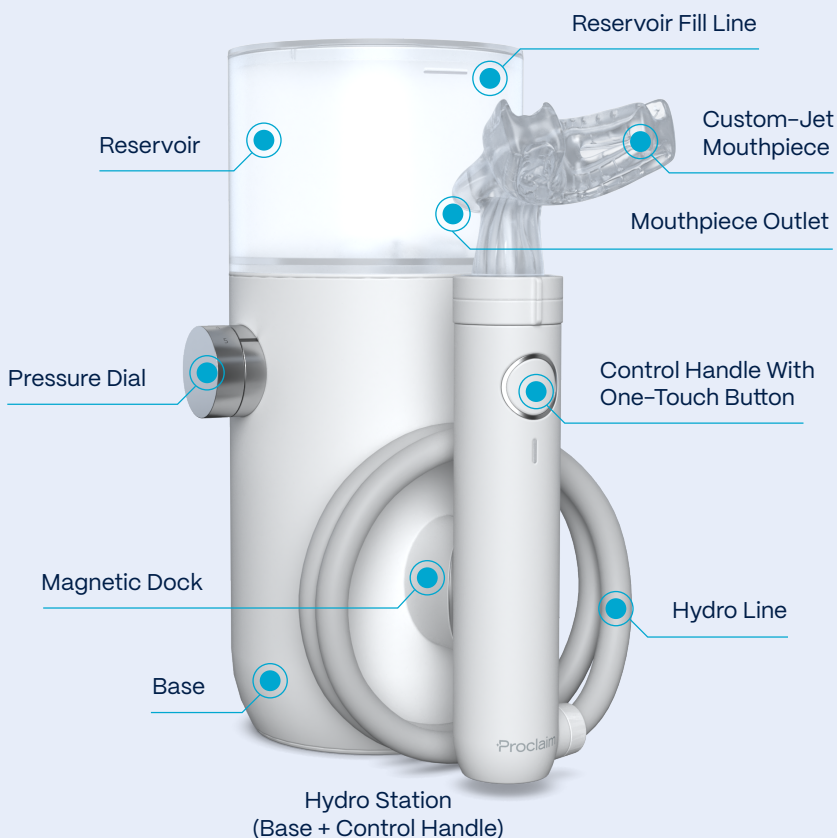


First Time Setup

- 1 Unpack the system and plug the Hydro Station into the outlet. Do NOT attach the Custom-Jet Mouthpiece yet.


NOTE: Handle your mouthpiece with care. It was created with customized technology and may be damaged if dropped. If dropped, inspect for damage. Always store the mouthpiece in its case to protect it.

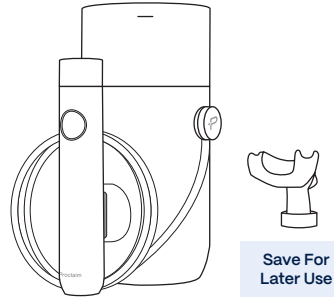

- 2 Remove the reservoir, fill it with **warm water** up to the fill line, and replace it on the base.

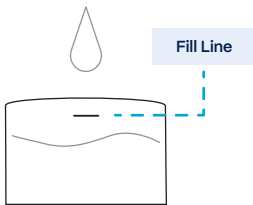


- 3 Set the pressure dial to **Level 1**.

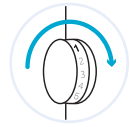


- 4 Remove the handle from the magnetic dock. Be careful not to accidentally push the button too early.

- 5 Remove the protective handle cap. You can save this and use it to cover the handle in between treatments.

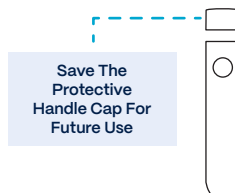


- 6 Point the end of the handle into the sink. Then press and hold the button until the cycle ends automatically.

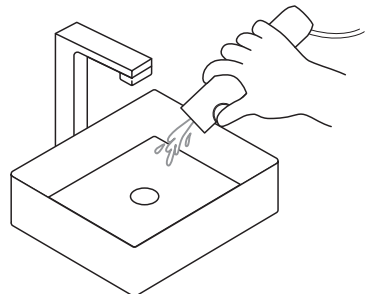

Prep For Treatment

Fill The Reservoir And Add Proclaim Oral Infusion


- 1 Fill the reservoir with **warm water** and place it back on the base. Make sure to fill it up to the fill line – the entire reservoir will be used for your treatment.

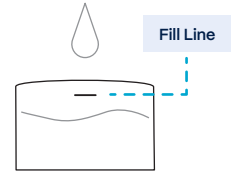

Reservoir Should Be Flush With The Base

- 2 For the most refreshing experience, add 20ml of the Proclaim Oral Infusion concentrate to the reservoir, using the included measuring cup.

NOTE: The Proclaim Oral Infusion has been specially developed for use in the Proclaim Custom-Jet Oral Health System. Other additives may cause damage to the device.

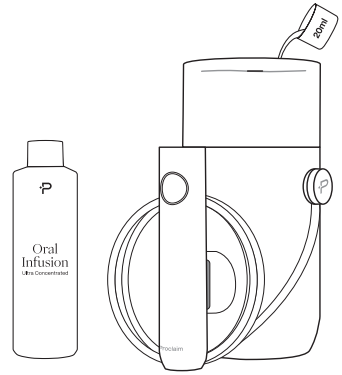


Elevate Your Proclaim Experience

For an invigorating minty fresh clean feeling with every use, don't forget to repurchase your Oral Infusion concentrate at proclaimhealth.com.



Set The Pressure Dial

Turn the dial to set the pressure level of your treatment.

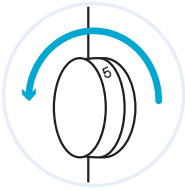


TIPS FOR A GREAT EXPERIENCE:


- ✦ Start at **Level 1** for your **first treatment**
- ✦ Increase over time to **Level 5**
- ✦ Continue all treatments at Level 5, the most effective pressure

Attach Your Custom-Jet Mouthpiece


- 1 Remove the handle from the magnetic dock.

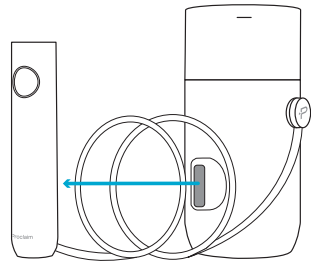

- 2 Attach your mouthpiece to the handle, positioning it so the bite surface is facing forward (aligned with the button).

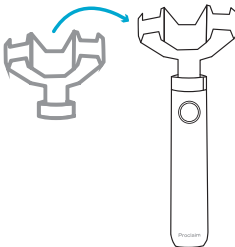

- 3 Twist the attachment ring clockwise until it clicks and the tab is centered.

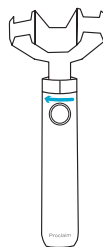

Attachment Ring
Should Click

Your 7 Second Treatment

Position Your Mouthpiece

- 1 Place your mouthpiece in your mouth and seal your lips around it.
- 2 Bite down securely without clenching—the mouthpiece material is hard and it may take a moment to adjust to the way it feels. The material enables precise targeting and effective water pressure.

TIP FOR A GREAT EXPERIENCE:

- ♦ Wet your custom mouthpiece with water before each use.



Run The Treatment Cycle

- 1 Lean over the sink. This way the liquid will drain out of the mouthpiece outlet into the sink.
- 2 Press and **hold the button on the control handle until the treatment is complete – it will stop automatically.** If you release the button early, the treatment will stop.



NOTE: At Level 5, the treatment will last about 7 seconds. At lower pressure settings, the treatment will last longer (up to about 25 seconds). The entire reservoir will be used in one treatment at all settings.

After Your Treatment Cycle

- 1 Detach your mouthpiece by twisting the attachment ring counterclockwise.
- 2 Rinse your mouthpiece under running water. Make sure to **store it in the provided case**, as it is precision technology and should be treated with care.
- 3 If you used the Oral Infusion concentrate, rinse the reservoir between uses.



For recommended cleaning and additional instructions, refer to the Proclaim Custom-Jet Oral Health System owner's manual at proclaimhealth.com/instructions or scan the code.

Still have questions? We are here to help.

Contact us at customercare@proclaimhealth.com



For additional information, instructions,
and FAQs, please visit our help center at
proclaimhealth.com/instructions
or scan the code.

